PL

EN

BEAUTIFLY

B-Glammi

P R O

series

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User manual

B-Glammi

P R O

series

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czujnik koloru skóry | 300 000 impulsów

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P R O

series

built-in hidden shaver	
skin color sensor	300 000 pulses
for different areas of the body (face, bikini, legs)	
5 levels of intensity	display
automatic pulse mode	

Hair removal technology using an IPL laser involves emitting a beam of light into the hair follicle and then destroying it. After exposing the bulb several times, it gradually breaks down.


This prevents hair regrowth and permanently removes it. B-Glammi has a built-in shaver, so you can do your scheduled IPL hair removal anytime, anywhere. In addition, the device has 5 adjustable levels of intensity, which means that you can adjust the device's operating power to the areas in which it is to be used.

Technical data

Rated voltage:	12V
Effective pulses:	300 000
Emitting frequency:	1-3.5 s/time
Emitted wavelengths:	550~1000nm
Optical exposure:	3.5J/cm ²
Optical areas:	3.22cm ²

Part description

1. Impulse window
2. Skin color sensor
3. On / off button and intensity button
4. Impulse button
5. Light indicator
6. Pulse intensity
7. Charging socket
8. Hair shaver
9. Charger
10. Goggles
11. Cleaning brush

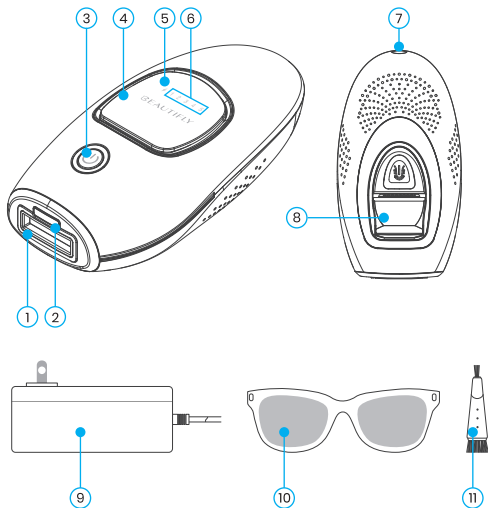Indications


No.	Item	Indication
1.	Shaving	Connect the device to power with adapter and the indicator light (white) flashes; gently push the groove upward with your thumb and the hair trimmer will open automatically to shave; return the groove back to stop shaving.
2.	Boot-up	After connecting the device to power, indicator light and all intensity lights flash once and then the indicator light continue to flash slowly. ↓ Press and hold on/off button for 1 sec. into mode 1, with the indicator light (white) and the corresponding intensity light (white) on.
3.	While working	When skin tone sensor detection is completed, the indicator light flashes, press flash button, the flash window flashes once; the indicator light flashes, and the corresponding intensity light (white) are on.
4.	Mode switching	While working, press on/off button (mode button) to switch to another mode. Emitting frequency is calculated according to the mode position, the indicator light (white) and the corresponding intensity light (white) are on. Mode 1, emitting frequency 0.8-1 s/time Mode 2, emitting frequency 1.5-2 s/time Mode 3, emitting frequency 2-2.5 s/time Mode 4, emitting frequency 2.5-3 s/time Mode 5, emitting frequency 3.5 s/time
5.	Power off	Press and hold on/off button for 1 sec. to stop working, the indicator light flashes slowly into standby; power off by unplugging the adapter.

Instruction of use


For the first use please follow the steps as below:

1. Wash and dry the skin where you will be using epilation. Do not use any creams and balms before performing the procedure. The skin must be dry and clean.

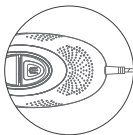2. Connect the device to the power supply.

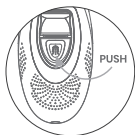3. Gently push the shaver switch upwards with your thumb and the razor will operate automatically.

4. Place the shaver at an angle of 90° on the skin and move it in the opposite direction to the hair growth.


5. Wear safety glasses attached to the device during laser hair removal.

6. Press and hold the on / off button for 1 sec to start the device. Place the flash window at an angle of 90°, pressing tightly against the skin.

7. Press the flash button to release a light pulse to remove hair. Then move the device sequentially place by place, releasing subsequent impulses.

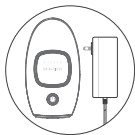8. Press and hold the on / off button for 1 sec. to turn off the device and disconnect the power supply.

9. Clean the device and the hair shaver (for more information, check the storage and cleaning sections).

Skin sensitivity test

If you are using the device for the first time or on a new area of the body, we recommend you take a skin sensitivity test.

Choose a **small** area near the area where you intend to do epilation. To do this, use the **lowest intensity level**.

Follow the instructions on how to perform the procedure below, but emit no more than 2-3 impulses. After the skin test, wait 60 minutes and check that the skin does not react.

If after 60 minutes the skin seems **normal** and there is no redness or other skin reaction, you are ready to perform the **full** procedure.

Before proceeding with the procedure

1. Before using B-Glammi, shave your hair from the skin. To do this, use the shaver located at the back of the device.
2. The shaver only works in standby mode. This means it will not work when you enable IPL.
3. The device has a skin tone identification system so that the light pulse is only released for the following skin tones. Typically, the larger the color difference between the skin and hair, the more effectively the hair is removed.

Skin tone available

4. When using, never fire the flash multiple times in the same area.
5. For large areas to be epilated, you can use the automatic pulse release mode by pressing and holding the pulse release button for 2 seconds. Then apply the device to the skin with the flash, press once and move place by place. The pulses will be triggered automatically.
6. If the skin will be exposed to the sun within 48 hours after use, apply SPF 50+ to the area of use.
7. After tanning, wait at least 48 hours before using the device.
8. The effect of hair removal varies from person to person due to different skin color and hair growth.
9. General schedule for hair removal during the full hair growth cycle:
 - From the first to the fourth week: twice a week;
 - From the fifth to the eighth week: once a week;
 - After eight weeks: once every half month or once a month.

[Warning]

1. The hair removal schedule varies depending on the type of hair. People with thick, thick hair should extend the epilation schedule.
2. Due to the high temperature of the device during laser hair removal, it is recommended not to disconnect the power supply immediately after continuous operation. Wait until the fan stops working (after 5 minutes) before disconnecting the power supply.

Contraindications

1. The device is not intended for children under 15 years of age. Young people aged 15-18 can use the device with permission and / or under the supervision of parents or guardians. Adults over 18 years of age can use the device freely.
2. Do not use if the skin in the area to be treated is tanned or burned. Using the device on tanned or burned skin can cause temporary or even permanent discoloration of the skin.
3. Do not use if you are pregnant, try to get pregnant, you suspect you may be pregnant or you are breastfeeding. The B-Glammi device has not been tested for this occasion, therefore, do not use it in these cases. Hormonal changes occurring during pregnancy or breastfeeding may increase the risk of skin damage or skin sensitivity.
4. Do not use the device if you have a disease such as epilepsy or or light-induced seizures or migraine. Energy flashes can cause severe headaches and migraines.
5. Do not use on areas covered with tattoos, dark brown or black places (such as dark freckles, birthmarks, or warts). Tattoos or dark spots can absorb too a lot of light energy that can cause swelling of the skin, temporary or permanent blisters, discoloration or whitening of the skin.
6. Do not use the device if you have a tendency to form keloids.
7. Do not use on the scalp, ears, neck, nipples or on mucous membranes (lips vulva smaller, around the vagina and anus). At those points of the bikini zone where the skin has stronger pigmentation, use a lower level of light intensity. Use in these places may cause swelling of the skin, blisters or temporary discoloration or skin bleaching.

8. Do not use if any active inflammation of the skin occurs in the treatment area, including acne, psoriasis, pimple, burns, herpes, abrasions or epidermis rash. These factors can make the skin more sensitive to light and make the procedure work unpleasant or cause swelling, blisters or infections on the skin. If the skin is enough healed, before using the device, consult a doctor.
9. Do not use the device if you are taking medicines that sensitize the skin to light including non-steroidal anti-inflammatory agents (e.g. aspirin, ibuprofen, cetaminophen), tetracyclines, phenothiazines, thiazides, diuretics, sulfonamides, DTICs, fluorouracil, vinblastine, griseofulvin, alpha-hydroxy acids (AHAS), Beta-Hydroxy Acids (BHAS), Retin-A, Accutane and / or local retinoids.
10. Do not use if you have a disease related to photosensitivity, e.g. porphyria, polymorphic light eruption, solar urticaria, lupus etc. These conditions make it the skin is more sensitive to light and the treatment can be unpleasant or causing swelling of the skin, blisters or infections.

Warning

-  Do NOT leave the device unattended while it is working.
-  Do NOT disassemble the device by yourself.
-  Do NOT use the device if the light tube is damaged.
-  Do NOT use or charge the device under the high-temperature and wet environment.
-  Do NOT use the product if it is hit or drops onto the ground heavily unless it is given a safety check.
-  Please ensure the power plug and cable are all in good condition in case of an electric shock.

* This list is not exhaustive. If there is any uncertainty as to the area of use or operation of the device, please consult the customer service department at the e-mail address: kontakt@beautifly.eu or by telephone on +48 12 345 09 72

Operating tips

1. Wear protective goggles when using the device.
2. Clean the area of skin to be treated before use. It is not allowed to apply any cosmetics to avoid the destroy of instrument and affect result.
3. Adjust the energy level from low to high level. On condition of safety, the higher the energy level is, the better result it will be. After treatment, it is normal to have skin redness, suggest to use ice pad after treatment.
4. After treatment then power off, please use wet cotton swab or cloth to clean treatment window. This can avoid hair and sebum of carbon sticking glass and affecting energy of shooting light.
5. Please don't let system and accessories touch liquid or put it in humid environment. Please store instrument in ventilated and dry environment and keep far away from children.
6. Please don't put metal object into heat emission hole. Otherwise there will be risk of electric shock.

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The crossed-out wheeled bin symbol indicates that in the European Union after the end of the product's useful life it must be disposed of in a separate, dedicated collection point. Do not dispose of these products together with unsorted general household waste

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Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beauty Factor Sp. z o. o. declares that the device equipment Beautifly B-Glammi is in compliance with Directive 2014/30/UE. The full text of the EU declaration of conformity is available at the following internet address:

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beauty Factor Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address biuro@beautifly.eu

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I represent that I have read and accepted the
terms and conditions specified in the Warranty
Certificate.

.....